



SkinTyte™

Deep heating for youthful skin appearance



 SCITON®


SkinTyte

the BEST WAY to TYTEN SKIN!

SkinTyte™ utilizes infrared energy to deeply heat dermal collagen while protecting the entire treated area with sapphire contact cooling. This process promotes contraction of the collagen, initiating the body's natural healing process. Following treatment, a new foundation of collagen is created leading to increased skin firmness.



BBL Smart Filters

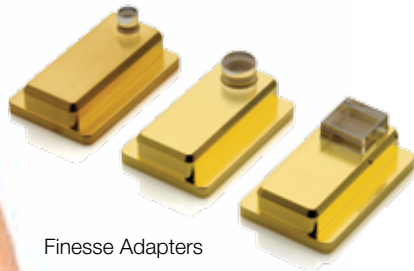

Finesse Adapters

SkinTyte For You:

SkinTyte delivers the fastest and most economic skin firming treatment solution on the market.

- **Excellent, Consistent & Predictable Results**

Uniform heating targets deeply and effectively to a desired area, immediately stimulating collagen remodeling and contraction.

- **Fast**

With a large 15 x 45 mm spot size and a rapid repetition rate, one can treat more patients in less time with the SkinTyte. For those smaller areas, magnetically attached Finesse Adapters provide added treatment flexibility.

- **No Consumables**

The SkinTyte has unlimited pulses under warranty and no consumables, which eliminates the extra cost of treatment tips.

“ Patients just love the SkinTyte because it causes no downtime and is painless and effective no matter where we've used it. It's been exciting to see the changes in the loose skin of the abdomen and underarms, as well as the face and neck. It's great to have no disposable cost and to cause no pain – unlike other tightening devices on the market.

Estelle Yamaki, M.D.

SkinTyte Results



Photos courtesy of Richard Sturm, M.D.



Photos courtesy of Gloria Lopez, M.D.

SkinTyte For Your Patients:

With the SkinTyte, you can offer your patients a highly requested, comfortable procedure with no downtime.

- **Powerful Cooling for Maximum Comfort**

Precision cooling with SkinTyte maximizes comfort and safety throughout the treatment. Only the BBL can cool the skin's surface to precise temperatures from 0°C to 30°C.

- **Safe for All Skin Types**

SkinTyte is safe for light and dark skinned patients who are interested in improving skin firmness.

- **No Downtime**

Patients are able to resume their regular activity immediately after with no visible signs of treatment.

Sciton's Post-Sale Support

As part of our commitment to your success, we provide not only the highest quality lasers, but also comprehensive training on your new laser, clinical application training, and practice marketing support programs. These include:

- **TopFlight™** (U.S. only) – we have partnered with experienced teams of third-party professionals to guide you through the integral aspects of building a cosmetic practice: including business development, marketing strategies, and product training and certification.
- **Practice Support Kit** (U.S. / International) – a comprehensive set of materials to assist you to effectively market applications to prospective patients as you seek to grow your practice.

Specifications

Wavelength	800 – 1400 nm
Pulse Width	1 sec – 15 sec
Speed	Up to 1.5 cm ² per sec
Skin Cooling	Adjustable from 0°C – 30°C
Spot Size	15 x 45 mm
Finesse Adapters	15 x 15 mm square, 11 mm & 7 mm round
Cooling Method	Continuous thermoelectric sapphire plates



visit us: www.sciton.com

925 Commercial Street, Palo Alto, California 94303
Phone: 888 646 6999 • Email: info@sciton.com
www.sciton.com

